Technical Data Sheet

PerCP-Cy[™]5.5 Mouse Anti-Mouse CD45.2

Product Information

Material Number: 552950

Alternate Name: Ly-5.2; T200; LCA; Leukocyte common antigen; Ptprc

 Size:
 0.1 mg

 Concentration:
 0.2 mg/ml

 Clone:
 104


Immunogen: B10.S mouse thymocytes and splenocytes


Storage Buffer: Aqueous buffered solution containing protein stabilizer and ≤0.09% sodium

azide.

Description

The 104 clone has been reported to react with CD45 (Leukocyte Common Antigen) on all leukocytes of most mouse strains (eg, A, AKR, BALB/c, CBA/Ca, CBA/J, C3H/He, C57BL, C57BR, C57L, C58, DBA/1, DBA/2, NZB, SWR, 129). This alloantigen was originally named Ly-5.1, and this was the designation at the time that the antibody was characterized. The designation was later changed from Ly-5.1 to Ly-5.2 to conform with the convention that the .2 alloantigen designations be assigned to the C57BL/6 strain. mAb 104 has been reported not to react with leukocytes of the mouse strains expressing the CD45.1 alloantigen (eg, RIII, SJL/J, STS/A, and DA). CD45 is a member of the Protein Tyrosine Phosphatase (PTP) family: its intracellular (COOH-terminal) region contains two PTP catalytic domains, and the extracellular region is highly variable due to alternative splicing of exons 4, 5, and 6 (designated A, B, and C, respectively), plus differing levels of glycosylation. The CD45 isoforms detected in the mouse are cell type-, maturation-, and activation state-specific. The CD45 isoforms play complex roles in T-cell and B-cell antigen receptor signal transduction. The 104 antibody has been reported to inhibit some responses of B cells, from mice expressing the CD45.2 alloantigen, to certain antigens and LPS. In addition, reduction of serum IgG levels and amelioration of autoimmune renal pathology were reported in mAb 104-treated systemic lupus erythematosus-prone mice.

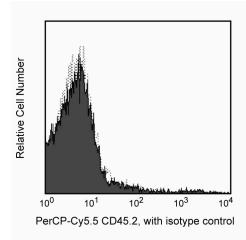

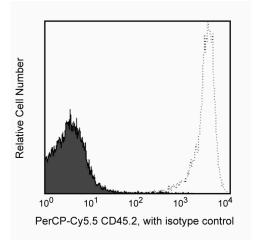

Differential expression of CD45.2 in SJL and BALB/c spleen. Splenocytes from SJL (left panel) and BALB/c AnN (right panel) mice were stained with either PerCP-Cy™5.5 Mouse IgG2a, κ Isotype Control (Cat. No. 550927; solid histograms) or PerCP-Cy™5.5 Mouse Anti-Mouse CD45.2 (Cat. No 552950/561096; open histograms). Flow cytometry was performed on a BD FACSCalibur™ flow cytometry system.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with PerCP-Cy5.5 under optimum conditions, and unconjugated antibody and free PerCP-Cy5.5 were removed. Storage of PerCP-Cy5.5 conjugates in unoptimized diluent is not recommended and may result in loss of signal intensity.

Application Notes

Application

Flow cytometry Routinely Tested

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 800.268.5430
 32.2.400.98.95
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



552950 Rev. 8 Page 1 of 2

Suggested Companion Products

Catalog Number	Name	Size	Clone	_
550927	PerCP-Cy [™] 5.5 Mouse IgG2a, κ Isotype Control	0.1 mg	G155-178	
554656	Stain Buffer (FBS)	500 mL	(none)	
553141	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block TM)	0.1 mg	2.4G2	
553142	Purified Rat Anti-Mouse CD16/CD32 (Mouse BD Fc Block TM)	0.5 mg	2.4G2	

Product Notices

- Since applications vary, each investigator should titrate the reagent to obtain optimal results.
- Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.
- PerCP-Cy5.5 is optimized for use with a single argon ion laser emitting 488-nm light. Because of the broad absorption spectrum of the tandem fluorochrome, extra care must be taken when using dual-laser cytometers, which may directly excite both PerCP and Cy5.5TM. We recommend the use of cross-beam compensation during data acquisition or software compensation during data analysis.
- PerCP-Cy5.5-labelled antibodies can be used with FITC- and R-PE-labelled reagents in single-laser flow cytometers with no significant spectral overlap of PerCP-Cy5.5, FITC, and R-PE fluorescence.
- Cy is a trademark of Amersham Biosciences Limited.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- Please observe the following precautions: Absorption of visible light can significantly alter the energy transfer occurring in any tandem fluorochrome conjugate; therefore, we recommend that special precautions be taken (such as wrapping vials, tubes, or racks in aluminum foil) to prevent exposure of conjugated reagents, including cells stained with those reagents, to room illumination.
- An isotype control should be used at the same concentration as the antibody of interest.

References

Greimers R, Trebak M, Moutschen M, Jacobs N, Boniver J. Improved four-color flow cytometry method using fluo-3 and triple immunofluorescence for analysis of intracellular calcium ion ([Ca2+]i) fluxes among mouse lymph node B- and T-lymphocyte subsets. Cytometry. 1996; 23(3):205-217. (Methodology: Flow cytometry) Johnson P, Greenbaum L, Bottomly K, Trowbridge IS. Identification of the alternatively spliced exons of murine CD45 (T200) required for reactivity with B220 and other T200-restricted antibodies. J Exp Med. 1989; 169(3):1179-1184. (Biology)

Morse HC 3rd, Shen FW, Hammerling U. Genetic nomenclature for loci controlling mouse lymphocyte antigens. Immunogenetics. 1987; 25(2):71-78. (Biology) Ogimoto M, Mizuno K, Tate G, et al. Regulation of lipopolysaccharide- and IL-4-induced immunoglobulin heavy chain gene activation: differential roles for CD45 and Lyb-2. Int Immunol. 1992; 4(6):651-659. (Biology)

Shapiro HM. Practical Flow Cytometry, 3rd Edition. New York: Wiley-Liss, Inc; 1995:280-281. (Methodology: Flow cytometry)

Shen FW. Monoclonal antibodies to mouse lymphocyte differentiation alloantigens. In: Hammerling GJ, Hammerling U, Kearney JF, ed. Monoclonal Antibodies and T-cell Hybridomas; Perspectives and Technical Advances. 1981:25-31. (Immunogen)

Shen FW, Tung JS, Boyse EA. Further definition of the Ly-5 system. Immunogenetics. 1986; 24(3):146-149. (Biology)

Yakura H, Ashida T, Kawabata I, Katagiri M. Alleviation of autoimmunity in BXSB mice by monoclonal alloantibody to Ly-5 (CD45). Eur J Immunol. 1989; 19(8):1505-1508. (Biology)

Yakura H, Kawabata I, Ashida T, Katagiri M. Differential regulation by Ly-5 and Lyb-2 of IgG production induced by lipopolysaccharide and B cell stimulatory factor-1 (IL-4). J Immunol. 1988; 141(3):875-880. (Biology)

Yakura H, Kawabata I, Shen FW, Katagiri M. Selective inhibition of lipopolysaccharide-induced polyclonal IgG response by monoclonal Ly-5 antibody. J Immunol. 1986: 136(8):2729-2733. (Biology)

Yakura H, Shen FW, Bourcet E, Boyse EA. On the function of Ly-5 in the regulation of antigen-driven B cell differentiation. Comparison and contrast with Lyb-2. J Exp Med. 1983; 157(4):1077-1088. (Biology)

BD Biosciences

bdbiosciences.com

United States
 Canada
 Europe
 Japan

 800.268.5430
 32.2.400.98.95
 0120.8555.90
 Asia Pacific Latin America/Caribbean 65.6861.0633

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violatio of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2014 BD



552950 Rev. 8 Page 2 of 2