Technical Data Sheet

BV421 Mouse Anti-Human CD38

Product Information

Material Number: 562445

Alternate Name: T10; ADP-ribosyl cyclase 1; Cyclic ADP-ribose hydrolase 1

 Size:
 25 test

 Vol. per Test:
 5 μl

 Clone:
 HIT2

Isotype:Mouse IgG1, κ Reactivity:QC Testing: Human


Workshop: III 155

Storage Buffer: Aqueous buffered solution containing BSA and ≤0.09% sodium azide.

Description

The HIT2 monoclonal antibody specifically binds to CD38. CD38 is a 45 kDa type II single-chain transmembrane glycoprotein present on thymocytes, activated T cells and terminally differentiated B cells (plasma cells). Other reactive cells include monocytes, macrophages, dendritic cells and some epithelial cells. The CD38 antigen acts as an ectoenzyme that catalyzes the synthesis and hydrolysis of a Ca++ mobilizing agent, cyclic ADP-ribose. This intracellular calcium plays an important role in cell signaling pathways. Reports describe CD38 as participating in adhesion with CD31, immunoregulatory functions involving signal transduction leading to cell growth, apoptosis, and differentiation.

The antibody was conjugated to BD Horizon™ BV421 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. With an Ex Max of 407-nm and Em Max at 421-nm, BD Horizon™ BV421 can be excited by the violet laser and detected in the standard Pacific Blue™ filter set (eg, 450/50-nm filter). BD Horizon™ BV421 conjugates are very bright, often exhibiting a 10 fold improvement in brightness compared to Pacific Blue™ conjugates.

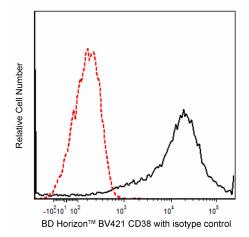

Flow cytometric analysis of CD38 expression on human peripheral blood lymphocytes. Human whole blood was stained with the BD Horizon™ BV421 Mouse Anti-Human CD38 antibody (Cat. No. 562444/562445; solid line histogram) or with BD Horizon™ BV421 Mouse IgG1, κ Isotype Control (Cat. No. 562438; dashed line histogram). The erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometry was performed using a BD FACSCanto™ II Flow Cytometer System.

Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV421 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV421 were removed.

Application Notes

Application

Flow cytometry Routinely Tested

Suggested Companion Products

Catalog Number	Name	Size	Clone	
554656	Stain Buffer (FBS)	500 ml	(none)	
555899	Lysing Buffer	100 ml	(none)	

BD Biosciences

bdbiosciences.com

United States Canada Europe Japan Asia Pacific Latin America/Caribbean 877.232.8995 800.979.9408 32.53.720.550 0120.8555.90 65.6861.0633 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only, Not for use in diagnostic or therapeutic procedures. Not for resale. Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



Product Notices

- 1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use 1 × 10⁶ cells in a 100-μl experimental sample (a test).
- 2. An isotype control should be used at the same concentration as the antibody of interest.
- 3. Brilliant VioletTM 421 is a trademark of Sirigen.
- 4. Pacific BlueTM is a trademark of Molecular Probes, Inc., Eugene, OR.
- 5. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
- Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
- For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at www.bdbiosciences.com/colors.
- 8. Please refer to www.bdbiosciences.com/pharmingen/protocols for technical protocols.

References

McMichael AJ, Beverly PCL, Gilks W, et al, ed. Leukocyte Typing III: White Cell Differentiation Antigens. New York: Oxford University Press; 1987. (Clone-specific)

Schlossman SF, Boumsell L, Gilks W, et al, ed. Leukocyte Typing V: White Cell Differentiation Antigens. New York: Oxford University Press; 1995. (Biology)

BD Biosciences

bdbiosciences.com

 United States
 Canada
 Europe
 Japan
 Asia Pacific
 Latin America/Caribbean

 877.232.8995
 800.979.9408
 32.53.720.550
 0120.8555.90
 65.6861.0633
 55.11.5185.9995

For country contact information, visit bdbiosciences.com/contact

Conditions: The information disclosed herein is not to be constructed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be help responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is stictly prohibited. For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



562445 Rev. 2 Page 2 of 2